



Tools needed for making wheelchairs

Ideas for setting up a workshop for *disabled* workers are discussed in Chapter 57 and p. 603 of Chapter 64. How you equip your workshop for making wheelchairs will depend on (1) how much money you have (or can borrow) to do it, (2) the kinds of chairs you hope to build (metal or wood), (3) the skills, *physical* and *mental* abilities, learning potential, and responsibility (regarding safety) of the workers, (4) the availability of electricity and power tools, (5) how many persons will be working, and (6) how many chairs you hope to produce.

Here we list the basic equipment you will need for making the 6 wheelchairs described in this chapter. Many choices are possible. More specialized parts of the work can be done by outside craftspersons. For example, in a wheelchair production center in Belize, axles must be machine tooled on a metal lathe. Local machine shops cooperate by doing this free.


CODE AN – Absolutely necessary N – A big help, but you might do without it (N) – Necessary only for axles ? – Depends on model	TYPE OF CHAIR					
	wood chair	re-bar a n d woven plastic	square metal tubes with wood seat and back	wheelchair with lying board	plywood	round metal tube
TOOLS REQUIRED	N	AN	N	AN	(N)	AN
bench vise	N	AN	N	AN	(N)	AN
tubing bender				AN		AN
welding (brazing) equipment	(N)	AN	N	AN	(N)	AN
metal saw	(N)	AN	AN	AN	(N)	AN
wood saw	AN			AN	AN	
hammer	AN	AN	AN	AN	AN	AN
wrench (set or adjustable)	N	N	AN	AN	N	AN
metal file and/or grinder	(N)	AN	AN	AN	(N)	AN
screwdriver	AN	AN	AN	AN	AN	AN
sewing equipment (hand or machine)			?	N?		N?
drill (hand or electric)	N	?	AN	AN	N	AN
drill bits for metal			AN	AN		AN
drill bits for wood	AN		AN		AN	
spoke wrench	?	?	N	N	?	N
bicycle pump	?	?	?	?	?	?
center punch	N	N	N	N	N	N
tape measure	N	N	N	N	N	N
carpenter's square	N	N	N	N	N	N

Terms for metal tube or bar used to build wheelchairs

- *Thin-wall* refers to thin steel tubing often used for electrical wiring work and sometimes for lightweight metal furniture.
- *Thick-wall* refers to heavy weight pipe such as the one used in plumbing.
- *Re-bar* refers to solid metal rod, usually used to reinforce cement.

Jigs or guides for more exact welding

For making the metal tube chairs and the welded wheel mounts and handrims of any of the chairs, your work will be easier and more exact if you make or purchase certain 'jigs' or guides to hold parts in the right place while you weld them. For example, to weld the front caster fork you can make a 'jig' like this. Details on 'jigs' and other techniques for making different wheelchair parts are well described in Ralf Hotchkiss's book *Independence Through Mobility* (see reference on p. 604). We strongly recommend it to any group planning to make wheelchairs.


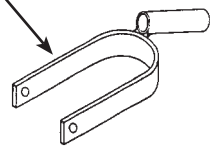


<b>Notes on measurements</b> For some of the wheelchair designs in this chapter, we give the measurements for a standard child's or adult's model. <b>Be sure to adapt the measurements to the size and needs of the particular child.</b>	In many countries inches (") are used for measurements of certain things, and centimeters (cm.) for others. We therefore also use both. <b>Centimeters</b> is abbreviated <b>cm.</b> and inches is abbreviated <b>"</b> . Two inches is written <b>2"</b> . <b>1"</b> equals 2.54 cm. You can use the scale on the edge of this page (and on the inside back cover) to change inches to cm.
---	---

# HEALTHLINK WOOD WHEELCHAIR

(Somewhat modified from Healthlink manual, see p. 604.)

The Healthlink wheelchair is built onto an ordinary **child's wood chair**. Measurements should be adjusted to the child's needs.

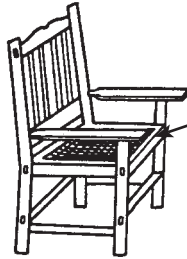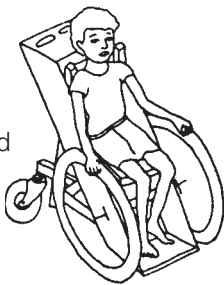


A webbed plastic seat lets air move through it and can be easily cleaned.

It uses standard 20" x 1 3/4" bicycle wheels and axles.



Basic carpentry tools are needed to build this wheelchair. It can be made in one day by someone with basic carpentry skills. The local blacksmith may be able to help weld together the wheel supports if you cannot. It is easy to add positioning aids or make other adaptations.

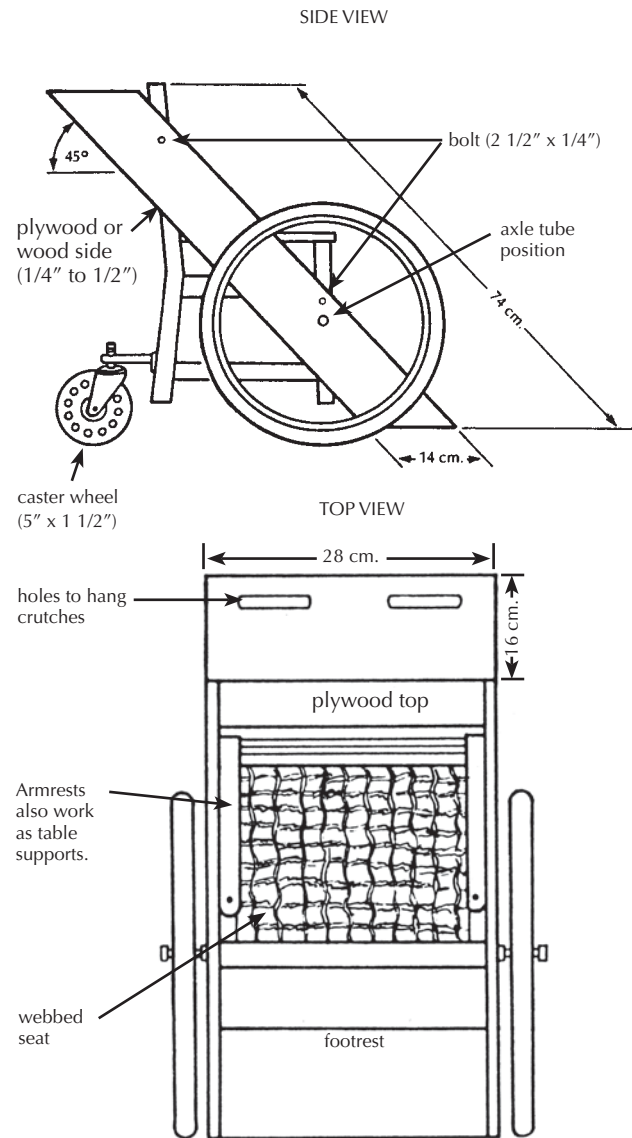

**AXLES**

Weld axles to ends of a steel tube 2 cm. longer than the chair is wide.

NO: Weld axle perfectly straight.

YES: Weld axle perfectly straight.


Pass axle tube through holes drilled through sideboards and front chair legs.

**FRONT VIEW**


**WARNING:** Use standard bicycle axles this way only for children under 20 kg (50 lbs.). A heavier child, or rough use, will bend or break the axle.

For children over 20 kg., use a stronger axle (see p. 623). Or support the bicycle axle from both sides (see p. 598).

## CASTERS

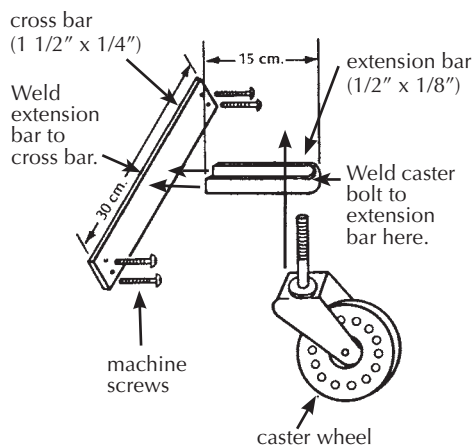


## BACK VIEW

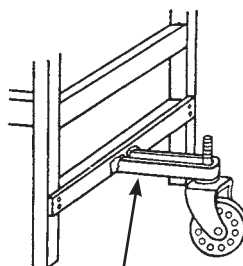


This backward extension distributes weight better and keeps chair from tipping backward on hills.

Two back casters are more stable and make it easier for the wheelchair to go up curbs.

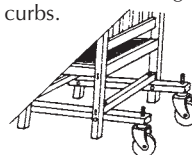

For **brake designs**, see p. 601 and 623. For other pictures and models of the Healthlink wheelchair, see p. 526, 592, 600, 601, 604, and 624.

## RE-BAR AND WOVEN PLASTIC WHEELCHAIR

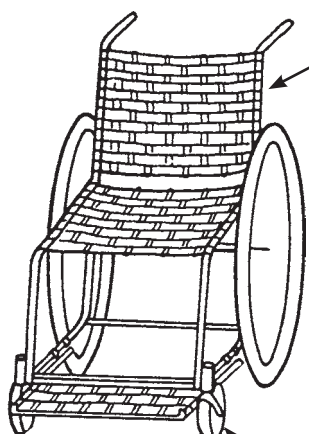


Weave back, seat, and footrest with ribbon, local wicker, cane rattan or polyethylene plastic.


Or use thin, lightly stretched strips of car inner tube, or canvas webbing.

For **front fork and casters**, use factory made casters or make your own (see above, and p. 597, 619, and 623).

## SIDE VIEW

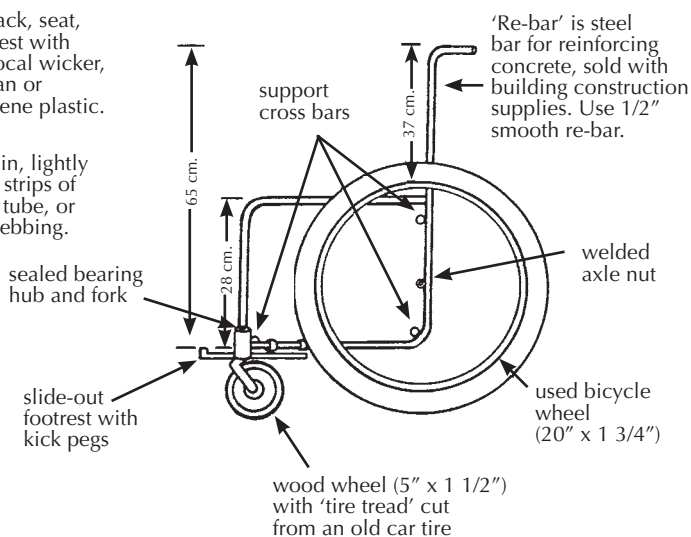

## SLIDE-OUT FOOTREST

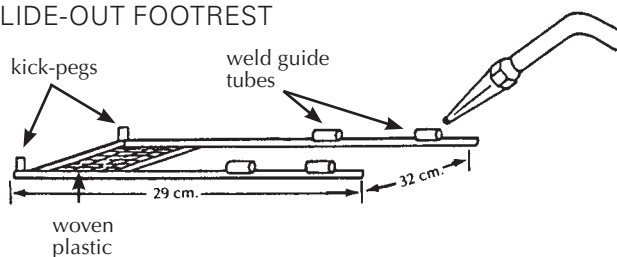


## MATERIALS NEEDED

- 1/2" re-bar (4 1/2 meters)
- inner tube strips
- bicycle wheels (2)
- front casters (2)
- webbing for seat


For **axle designs** see p. 597, 598, 615, and 623.

## SQUARE TUBE WHEELCHAIR

This wheelchair, like other steel tube chairs, should use only thin-wall tubing. To keep costs down, check with various sources of materials and ask at small fix-it shops for advice and possibly even some free scrap material. Metal scrap heaps are great for materials.

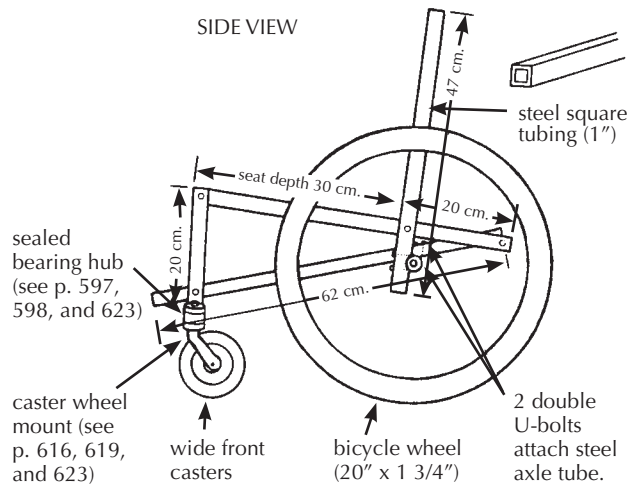

FRONT VIEW

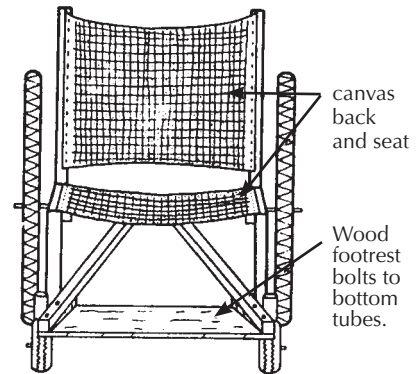


TOP VIEW

Side U-bolts secure axle tube.

A weld here adds strength, but is not necessary.

Threaded rod connects bottom tubes to back tubes.

Middle U-bolts secure axle tube.

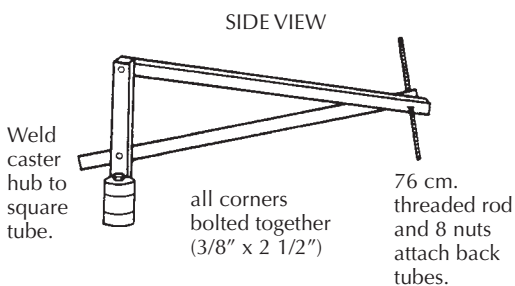


### MATERIALS NEEDED

- thin-wall square tubing (1" x 3.64 meters)
- thick canvas cloth (1 square meter)
- galvanized steel tube (1/2" x 66 cm.)
- bicycle wheels (2) (20" x 1.75")
- caster wheels (2) (wood or rubber)
- threaded rod (3/8" x 38") (Use extra 20" to bend 4 U-bolts.)
- 2 front casters
- 21 3/8" nuts and 12 screws for seat and back supports

### HOW TO MAKE YOUR CHAIR

1. Review drawings. Adjust measurements to fit child.
2. Cut all sections of square tubing. Make sure that matching tubes are equal in length.
3. Drill holes in bottom tubes and pass the threaded rod through them. Adjust nuts until a 'V' is formed. (Weld tip of 'V' for extra strength.)
4. Drill all holes in seat tubes. Pass threaded bolt through seat holes.
5. Drill holes in back support tubes and front caster tubes. Bolt to frame.
6. Weld axle nuts to ends of axle tube. Drill holes for U-bolts and bolt axle tube to frame.
7. Weld front caster forks to front tubes.
8. Sew cloth back and seat supports. Screw into place.
9. Cut out and bolt wood footrest to frame. (Use wedges to get the angle right.)
10. Attach axle tube with U-bolts and put on the wheels.
11. Paint frame to help keep tubes from rusting (if not galvanized).

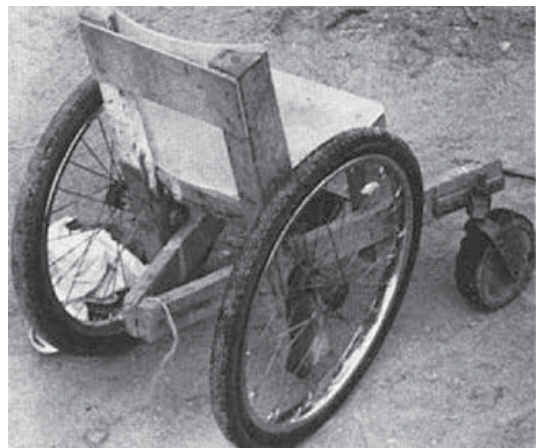


The same design can be made of wood.



WHEELCHAIR WITH LYING BOARD


This is useful for an active child who must lie face down to heal pressure sores or to stretch hip and knee contractures.

The board is sloped so that the child can play, look ahead, and move about more easily. If necessary, you can make the lying board adjustable so that the child can rest lying flat. This helps to improve *circulation* and to prevent swelling of the feet.

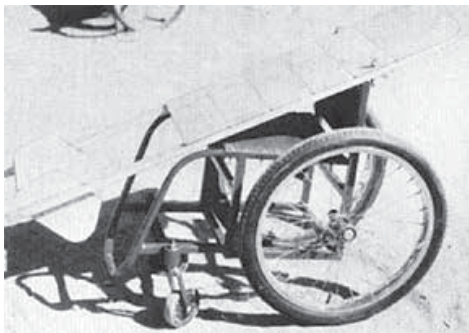


After the pressure sores heal, the lying board can be removed and the frame is easily adapted to form a lightweight wheelchair.

The design we show uses a simple, non-folding steel tube wheelchair frame with a wooden lying board mounted on top. However, many other designs are possible. (See, for example, the photo of a lying and standing wood wheelchair on p. 190.)

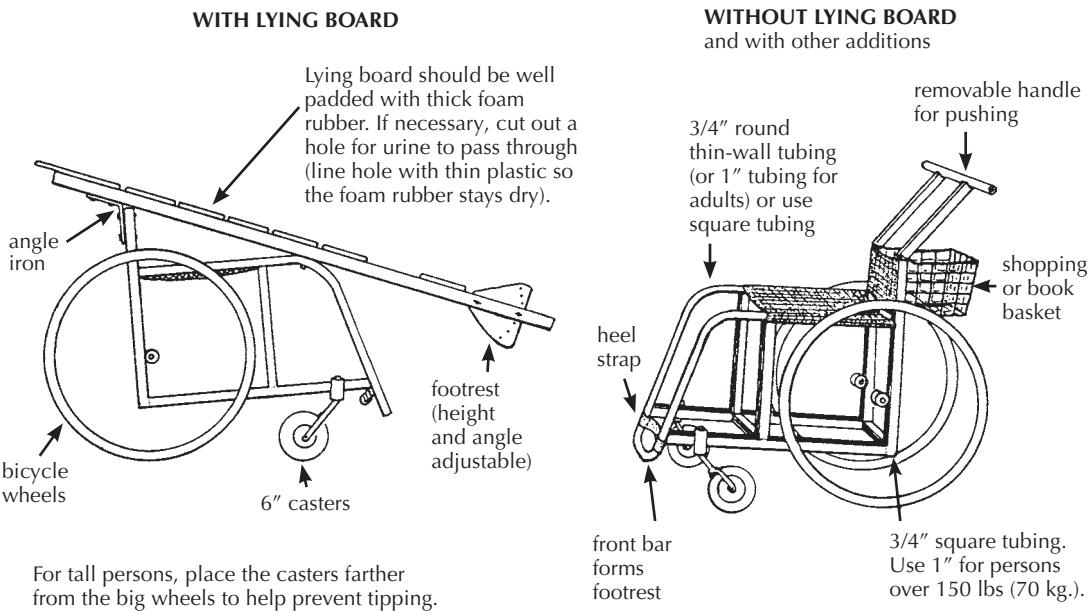


THE LYING BOARD

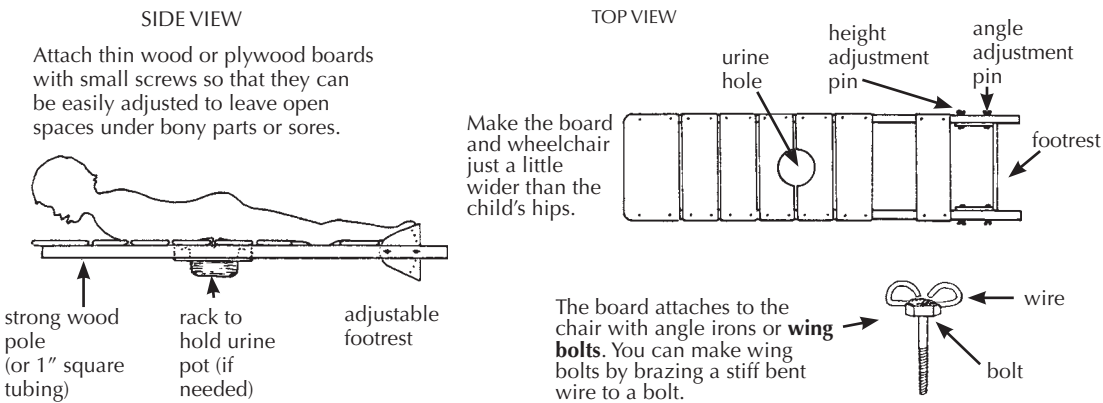


FOOTREST


Use thin wood or plywood. (Pad sides and bottom well to prevent sores. Examine feet daily.)

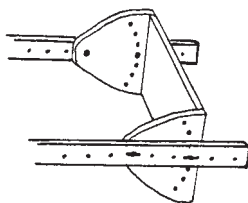

REMOVABLE HANDLE

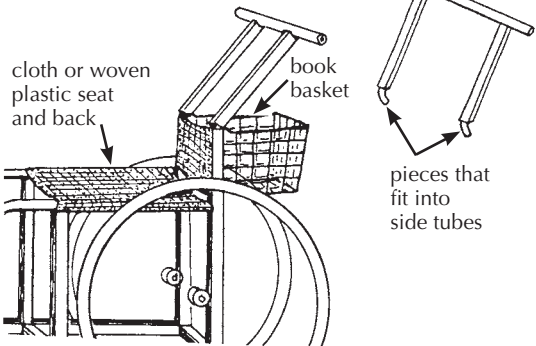

FRONT CASTER WHEEL

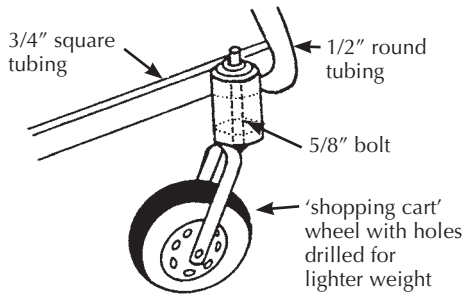



You should now have enough information to make a wheelchair with a lying board without step-by-step instructions. Adapt it, and make it the size to fit the child that needs it.



Wheelchair with lying board. A wide strap holds the child in place (but take care it does not press on sores).



Wheelchair without lying board.

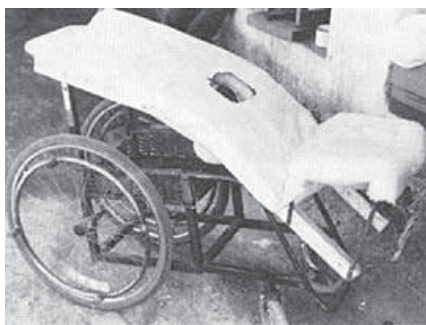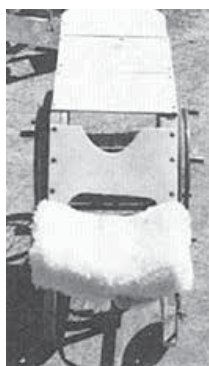

A variation of the wheelchair with lying board (p. 618) adapted for a paraplegic child with both contractures and pressure sores of his hips and knees. Urine is collected in a plastic container. The wheelchair seat has been converted into a basket.

**CAUTION:** Remember that a child who has some pressure sores can easily get new ones. Be sure the child lies and sits so that there is little or no pressure over bony places. **Examine her whole body at least once a day and try to keep her dry.**

CP **PLYWOOD FRAME WHEELCHAIR**


This can be easily built by someone with basic carpentry and welding skills. Positioning aids (head rest, hip pads, etc.) can be easily added. The chair can be designed to meet a child's particular needs. For example, if the child sits well without extra support, the tops of the side pieces can be removed to allow more freedom of movement.

A plywood frame is a low-cost alternative to metal. However, if not made well, or if left out in the rain, the chair may weaken and the plywood can split. As with any wheelchair, it must be protected from misuse, periodically examined for weaknesses, and promptly repaired.

For active children the wheelchair can be strengthened by reinforcing all joints and by adding strong hubs and axles (see p. 623).



See model on p. 621 →

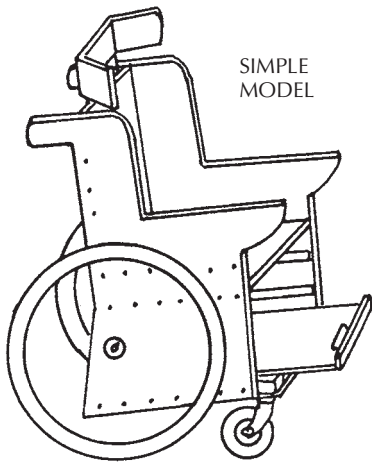


**HOW TO MAKE YOUR CHAIR**

1. Review drawings of chair and adaptive equipment.
2. Cut out the two side pieces to the same shape; sand with sandpaper.
3. Cut out back support, seat, and bottom piece of chair; sand with sandpaper.
4. Screw or nail seat and bottom piece to back piece.
5. Screw or nail side pieces to seat, bottom, and back.
6. Check that all pieces are lined up straight. Then add glue and more screws or nails for strength.
7. Cut out footrest and guide brackets for footrest.
8. Screw or nail guide brackets to side pieces under seat.
9. Bolt front casters to chair and assemble rear axle tube.
10. Drill holes in side pieces for axle tube; mount tube and rear wheel.
11. Let glue dry 1 to 2 days; check for strength of all wood joints.

These measurements are for a 4 to 8-year-old child.

- MATERIALS NEEDED**
- 3/8" plywood (1 sheet)
  - 20" bicycle wheels (2)
  - small caster wheels (2)
  - 1/2 steel tube (66 cm. long)
  - wood glue
  - sandpaper
  - screws
  - nails
  - 1/2" by 1/4" wood strips (6 x 46 cm. long)

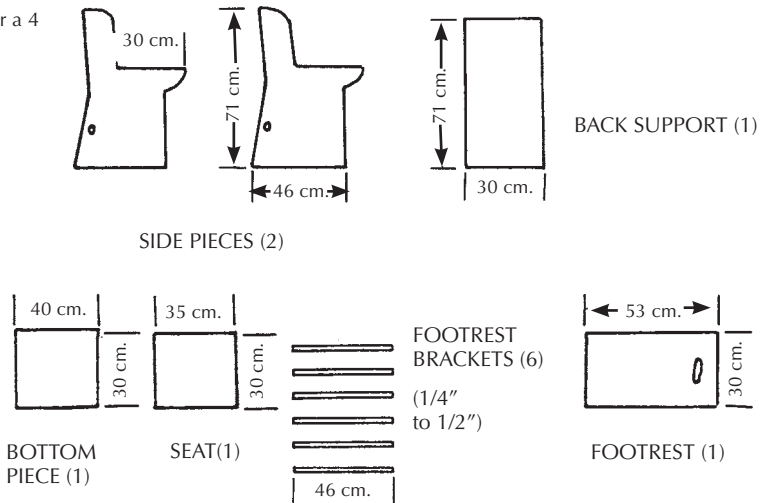






## A plywood wheelchair with many adaptations

This wheelchair has a variety of additions sometimes needed for a small child who has poor body control, head control, and urine or bowel control. The head support and armrests fit into wooden holders and can be easily removed. A lap table can be easily added. Holes can be cut out for chest and hip straps for extra support.

SIDE VIEW

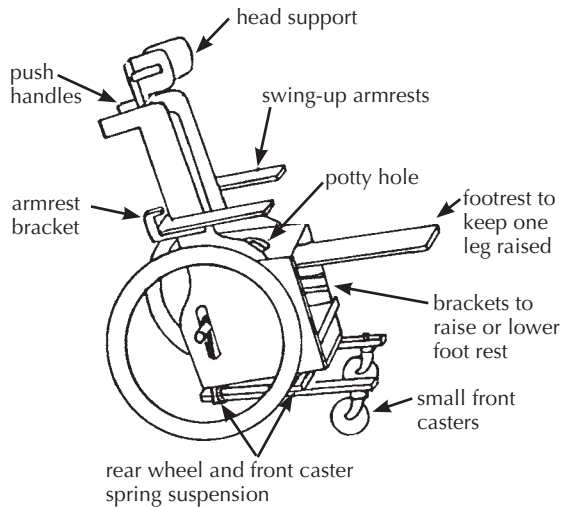

BACK VIEW

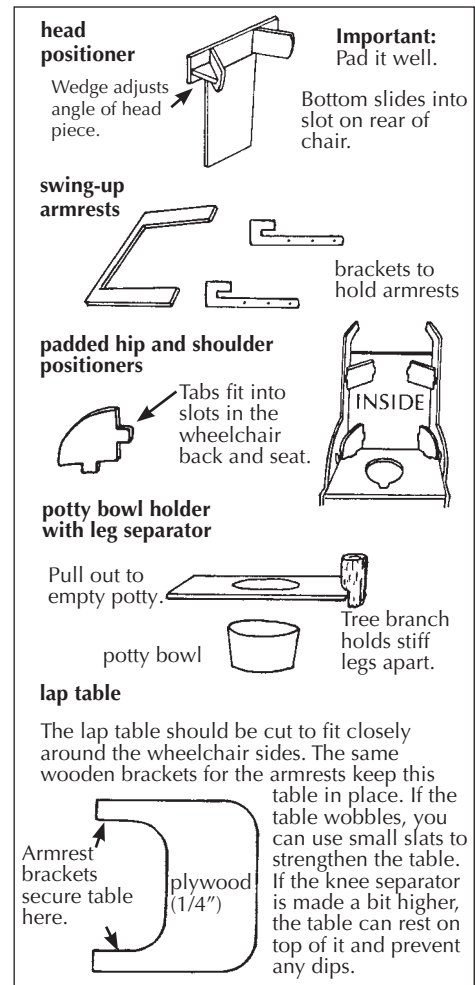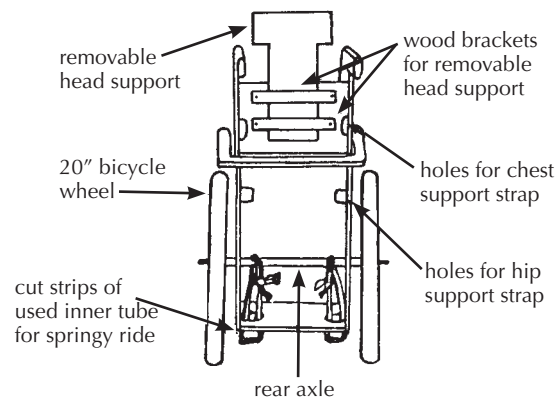


## SPRINGS FOR ALL 4 WHEELS

This plywood wheelchair has a springy ride. Old inner tube rubber strips connect the rear wheel axle to the wood strips holding the front caster wheels. These wooden strips should be strong enough to withstand the springy motion of the front casters.

Special cut-away slots allow the rear axle to move up and down freely. Other cut-away slots in the bottom of the wheelchair allow for the inner tube strips to be wrapped around the wooden caster strips. The tighter the inner tube strips are wrapped, the less bouncy the ride becomes.

To build your own strong rear hub and axle, see p. 623. If you want to use hubs from bicycle wheels, see p. 597.

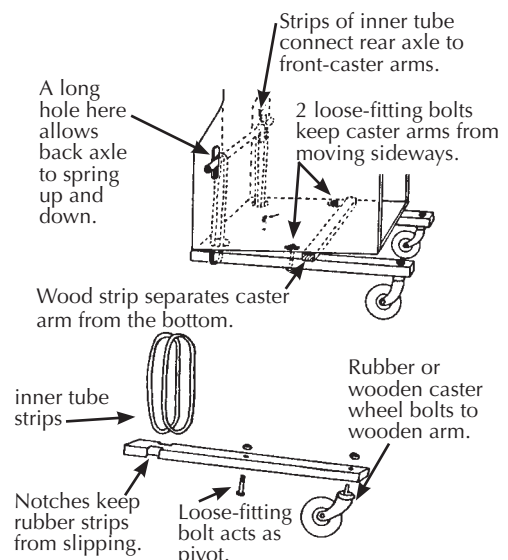

WHIRLWIND STEEL TUBE WHEELCHAIR



The whirlwind (ATI-Hotchkiss) wheelchair is a very strong lightweight folding chair. On rough ground it rides more easily and lasts longer than more costly factory-made chairs. If it breaks, it can be fixed by the neighborhood metalworker. It is narrow and helps the rider to move about crowded rooms.

The frame of this chair is made of thin-wall steel tubing that is easy to shape by someone with basic mechanical and welding skills. It can be built in about 4 days in a small metalworking shop. More than 10 groups of disabled mechanics throughout Latin America are building this wheelchair—often at less than a quarter the cost of imported wheelchairs.

Most materials for this chair can be obtained locally. It uses standard 24" (or 26") bicycle wheels. The extra strong hubs (see p. 623) use standard small machinery bearings (which can often be obtained used for free or at low cost from electric machinery repair shops). The axles are 5/8" (1.6 cm.) steel bolts. Seating is canvas (heavy cloth). If the small front wheels are not available, you can make them out of wood (see p. 597 and 616).

The curved fender bar that follows the shape of the tire makes transfers easier. The lightweight folding footrests are narrow at the front, for moving more easily in crowded spaces.

Plans for making hubs, casters, and brakes are on the next page. Complete plans for making this wheelchair are in the book *Independence Through Mobility* (see p. 604). The book is essential for anyone planning to build this chair.

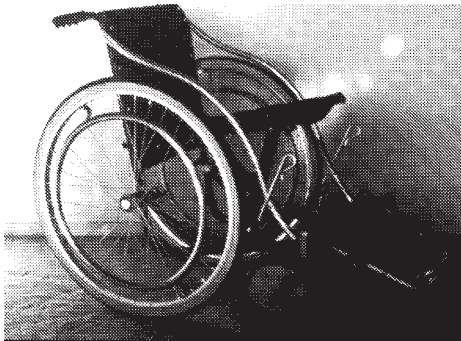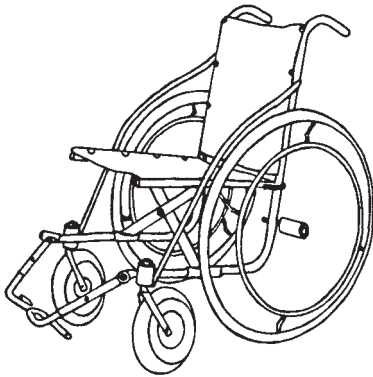


Model with wooden front wheels

MATERIALS NEEDED



- thin-wall tubing (from 1/2" to 1 1/4")
- thick-wall tubing (5/8" inside diameter)
- thick canvas or nylon cloth (2 meters)
- square tubing (thin-wall)
- bicycle rims and spokes (24" or 26" diameter)
- caster wheels (2)
- used sealed bearings (8)
- re-bar steel (3/8" round)

- flat bar steel (1/16" x 3/8")
- axle bolts (4) (5/8" x 5")
- washers (4) (1" diameter, 16 upholstery)
- screws (8 upholstery)
- machine screws (8) (1/4" x 1 1/2")
- paint or chroming chemicals
- bronze welding rod, flux
- bicycle tires and inner tubes (24")

FOLDING FOOTREST

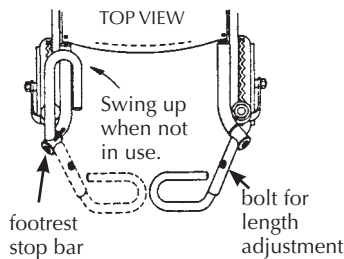

X-BRACE

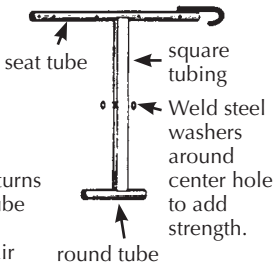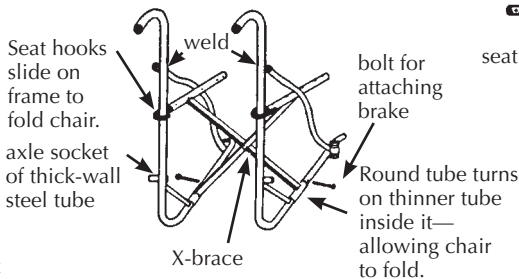


For a photo of this chair, see p. 536.


## DETAILS OF HOW TO MAKE WHEELCHAIR PARTS

(can be used with many wheelchair designs)


### HEAVY-DUTY WHEEL HUBS

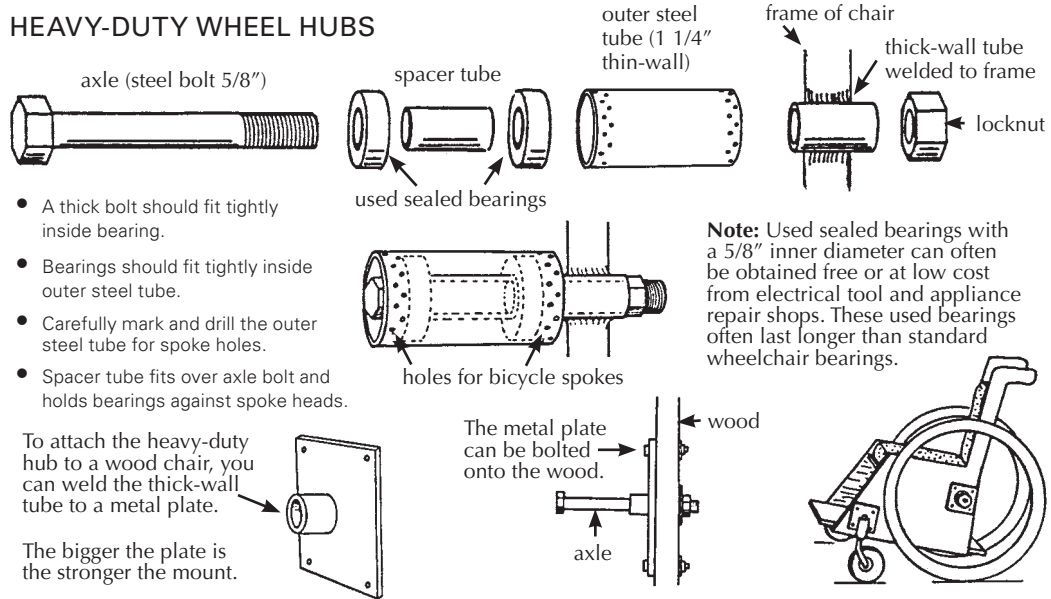


### CASTER FORK AND HUBS

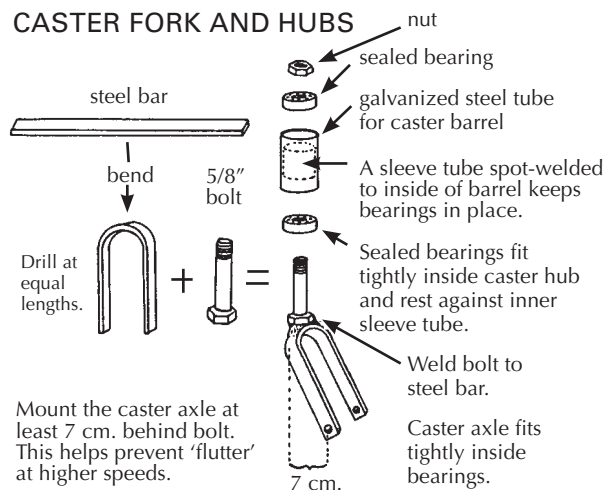

### BRAKES

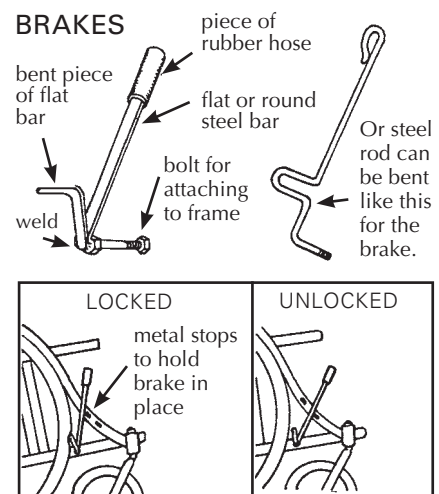

### SEAT AND BACK MADE OF CANVAS (STRONG CLOTH)

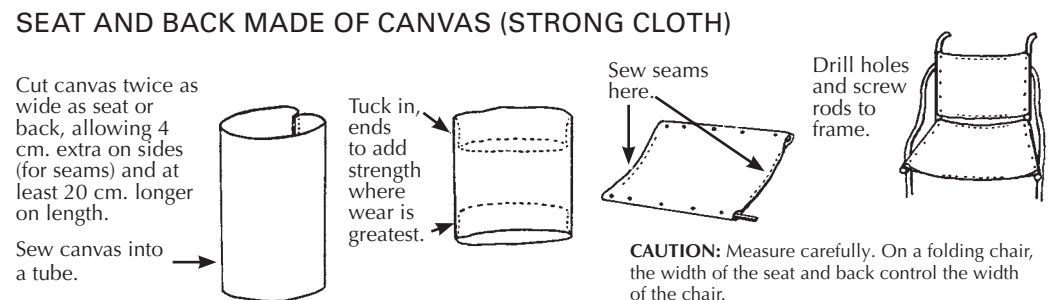

**For designs of other wheelchair parts, see the following pages:**

wheels: 594, 596, 597, 616, 619  
seats and backs: 595, 615, 616, 617, 619, 620  
tires: 596  
armrests: 599, 621

footrests: 600, 616, 619, 621, 622 axle  
mounts: 597, 598, 615  
handrims: 601 cushions:  
200, 609



Examples of locally made wheelchairs


CP




The plywood wheelchair on p. 620, with the armrest in place (left) and swung back (right).




A plywood wheelchair for a child with cerebral palsy with inner tube stretching aids to gently pull his feet and straighten his severe knee contractures.

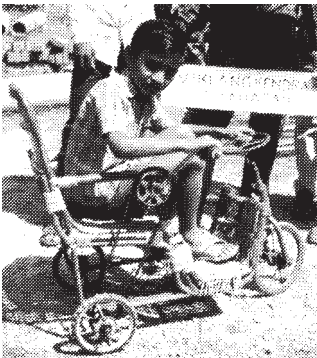


A bamboo hand-powered tricycle made at Viklang Kendra (People's Village), Allahabad, India.

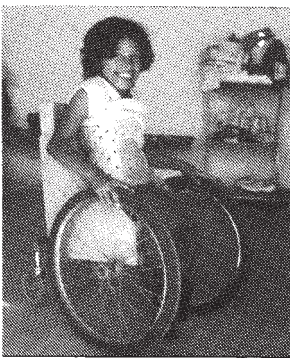


A wheelchair made completely of paper, including the wheels. Paper is glued together using rice flour in water (Zimbabwe).

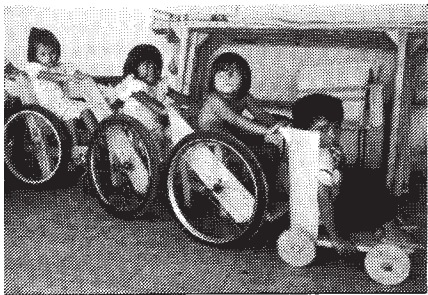


A wood design of the wheelchair on p. 617, two Healthlink wheelchairs, and a 'trolley' made from half of a plastic bucket and wood wheels.

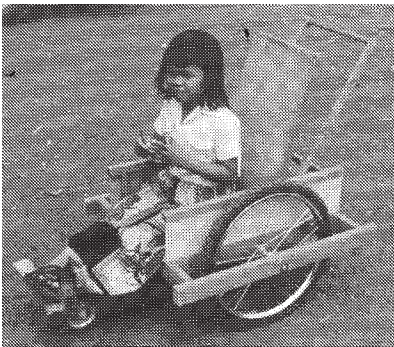

A wood wheelchair in Thailand. The bicycle wheel axles are supported on both sides to keep them from bending.

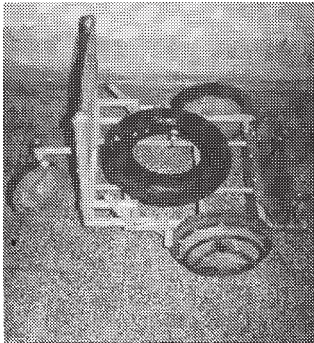

A metal frame, wood wheel 'trolley' in Bangladesh (see p. 572). The rubber tube serves as a cushion and also as a toilet seat.

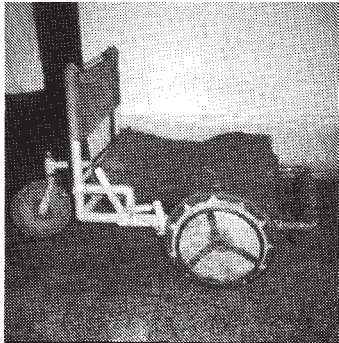

This trolley, also from Bangladesh, uses a cushion made of coconut fiber covered with rubber (see p. 199).

For more examples of wheelchair designs, see p. 65, 86, 98, 189, 190, 229, 288, 343, 430, 441, and 526.